


References:

1. Loveday HP, Wilson JA, Pratt RJ et al. (2013) Epic3: National Evidence-Based guidelines for Preventing Healthcare-Associated Infections in NHS Hospitals in England. *Journal of Hospital Infection*; 86S1: S1-S70.2
2. National Institute for Health and Care Excellence (2012) Guidance 139. *Infection: Prevention and Control of Healthcare-associated Infections in Primary and Community Care*. London: National Institute for Health and Care Excellence; available at www.nice.org.uk
3. Healthfirsteurope.eu. 2021. INCREASING ADHERENCE TO CAUTI GUIDELINES: RECOMMENDATIONS FROM EXISTING EVIDENCE. [online] Available at: <http://healthfirsteurope.eu/wp-content/uploads/2021/06/Increasing-adherence-to-CAUTI-guidelines-recommendations-from-existing-evidence.pdf> [Accessed 10 August 2021]

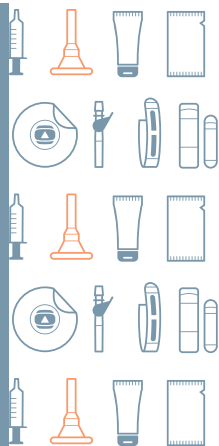

Optimum Medical Solutions Ltd
Tennant Hall, Blenheim Grove,
Leeds, LS2 9ET, United Kingdom
T: +44(0) 845 643 5479
enquiries@optimummedical.co.uk
www.optimummedical.co.uk

EC REP

MT Promed Consulting GmbH
Ernst-Heckel-Straße 7
66386 St. Ingbert Germany

UROLOGY RANGE



optimummedical®

| constantly creating better

We understand the importance of having the right product. We know that for patients who need a catheter, independence and confidence are key. As specialists in healthcare lubrication and urology, we've spent many years listening to clinicians to truly understand the needs of both them and their patients. Then doing what we do best - constantly creating better. Better products that offer quality, care and comfort to everyone that needs them.








Urology - Lubrication



optimummedical®

Contents

Procedure Packs	5
Ugo Pro (catheterisation procedure pack)	6
The Ugo Box	8
Drainage	11
Ugo Foley Catheter Kit	12
Ugo Leg Bags	14
Ugo Pro Leg Bags	16
Ugo 2L Drainage Bags	18
Ugo Sheath	20
Ugo Catheter Valve	22
Ugo Stand with Dignity	24
Fixation Devices	26
Ugo Fix Gentle	28
Ugo Fix Catheter Strap	30
Ugo Fix Sleeve (leg bag holder)	32
Ugo Fix Bag Belt	34
Multi Fix Gentle	35
Catheter Drainage and Management Pack	36
Ugo 4 Weeks	38
Ugo Home	40
Catheterisation Jelly	42
OptiLube Syringes	44
OptiLube Active	45
OptiLube Active CHG Free	46
Optilube	47
Choosing the right catheterisation jelly for your patient	48
Sustainability	50

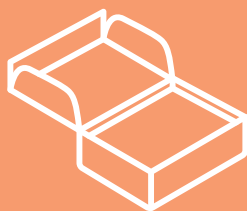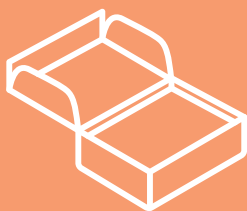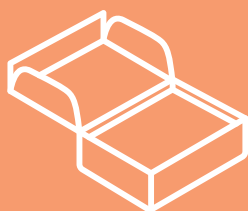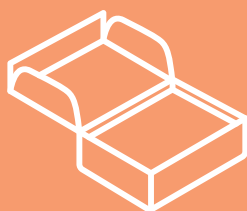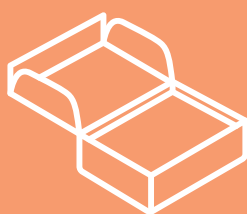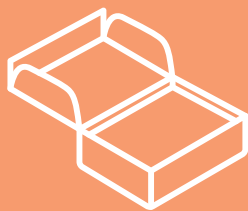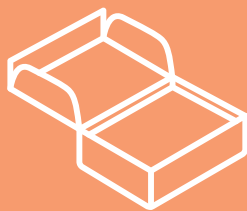

Procedure Packs

Our procedure pack range has been designed in line with EPIC¹ guidelines and clinician input to support best practice aseptic catheterisation, promote efficiency, and provide a cost-effective solution in one handy pack.

A study by the European Association of Urology Nurses (EAUN)³ found an average reduction of around 80% in CAUTI when urinary catheterisation kits/sets were introduced in a primary care setting.

IN THIS RANGE

Ugo Pro	6
The Ugo Box	8



Developed with clinician input, the **Ugo Pro** provides everything needed to support best practice insertion and removal of transurethral or suprapubic catheters in a community setting. Each pack contains a pre-wrapped sterile field and the essential routine components required to carry out the procedure, helping to minimise the risk of infection.

The **Ugo Pro** helps to minimise the unnecessary waste which would accumulate if each item were to be obtained separately and assists in streamlining stock management.

UGO PRO CONTAINS:

- 1 Clear instructions for use.
- 2 1 x outer wrapper – creates an aseptic field for the components to be placed upon.
- 3 1 x disposable apron.
- 4 1 x waste bag.
- 5 1 x grey waterproof drape – protects the bed from spilt fluids during procedure. Assists with aseptic technique.
- 6 1 x blue waterproof “O” drape - provides dignity during the procedure. Assists with aseptic technique.
- 7 1 x Ugo 10ml syringe – to remove fluid from Foley catheter balloon prior to catheter removal.
- 8 6 x non-woven swabs – to assist with removal and general catheterisation procedure technique.
- 9 5 x non-woven balls – for cleansing catheter insertion site during procedure (once coated in saline).
- 10 1 x 20ml sterile saline (0.9% w/v sodium chloride) - to empty over non-woven balls for cleansing catheter insertion site.
- 11 1 x clear tray with 2 compartments.
- 12 3 x pairs of medium nitrile **Ugo Pro** gloves – assists with best practice technique.

Ugo Pro (catheterisation procedure pack)


PR CODE	PR NAME	OUTER CASE
3400	Ugo Pro (Catheterisation procedure pack)	Case of 80



Key features:

- The unique design of the **Ugo Pro** allows freedom to choose which size/type of catheter, leg bag or catheter valve, and lubricant to use depending on patient or clinician preference.
- **Ugo Pro** includes sterile saline, which is recommended by EPIC3 as cleaning solution.
- **Ugo Pro** comes with 3 pairs of gloves, which is the required amount for best practice catheter removal and insertion.
- The drapes included in the **Ugo Pro** are different colours for easy differentiation.



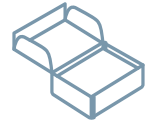

The **Ugo Box** is a fully equipped catheterisation kit, containing everything needed to safely remove and insert a transurethral or suprapubic catheter in a community setting, following best practice catheterisation technique.

The **Ugo Box** gives you, the clinician, the option to choose the catheter size (CH), catheterisation jelly syringe and leg bag, ensuring you use the most suitable catheterisation kit for your patient. Available in 12 variations.

THE UGO BOX CONTAINS:

- 1** 1x **Ugo Pro** (catheterisation procedure pack with 20ml sterile saline and 10ml syringe)
- 2** 1x **Ugo Foley Catheter** (all silicone)
- 3** 1x **OptiPure** (10ml sterile water syringe)
- 4** 1x 11ml catheterisation jelly syringe (**OptiLube** or **OptiLube Active CHG Free**)
- 5** 1 x **Ugo Leg Bag**
Choose either:
Ugo Leg Bag 1C - short tube, 500ml, lever tap.
Ugo Leg Bag 2C - long tube, 500ml, lever tap.
- 6** 1x **Ugo Fix Gentle** (catheter clip)
- 7** 1x pair of **Ugo Fix Leg Bag Straps**

The Ugo Box




PR CODE	UGO FOLEY CATHETER CH SIZE	UGO LEG BAG	11ML CATHETERISATION JELLY SYRINGE	SUITABLE FOR:	OUTER CASE
3450	12CH	Ugo 1C	Optilube Active CHG Free	Transurethral catheterisation only.	Case of 9
3451	14CH	Ugo 1C	Optilube Active CHG Free	Transurethral catheterisation only.	Case of 9
3452	16CH	Ugo 1C	Optilube Active CHG Free	Transurethral catheterisation only.	Case of 9
3453	12CH	Ugo 1C	OptiLube	Transurethral or suprapubic catheterisation.	Case of 9
3454	14CH	Ugo 1C	OptiLube	Transurethral or suprapubic catheterisation.	Case of 9
3455	16CH	Ugo 1C	OptiLube	Transurethral or suprapubic catheterisation.	Case of 9
3456	12CH	Ugo 2C	Optilube Active CHG Free	Transurethral catheterisation only.	Case of 9
3457	14CH	Ugo 2C	Optilube Active CHG Free	Transurethral catheterisation only.	Case of 9
3458	16CH	Ugo 2C	Optilube Active CHG Free	Transurethral catheterisation only.	Case of 9
3459	12CH	Ugo 2C	OptiLube	Transurethral or suprapubic catheterisation.	Case of 9
3460	14CH	Ugo 2C	OptiLube	Transurethral or suprapubic catheterisation.	Case of 9
3461	16CH	Ugo 2C	OptiLube	Transurethral or suprapubic catheterisation.	Case of 9

Please note: The **Ugo Foley Catheter** balloon size is 5 - 10ml.
OptiLube Active CHG Free is not suitable for use during suprapubic catheterisation procedures.
Ugo Foley Catheter 12CH features a non-integrated balloon.
Ugo Foley Catheter 14CH and 16CH feature an integrated balloon.



Key features:

- Suitable for the insertion and removal of transurethral and suprapubic catheters.
- Recyclable packaging.
- Reduces waste - by providing all the items required in one box, this reduces stockpiling, over-prescribing and excessive packaging.
- Promotes efficiency - reduces time spent setting up the procedure, resulting in more time to spend with the patient.
- Latex-free.

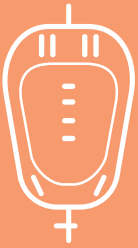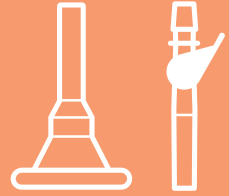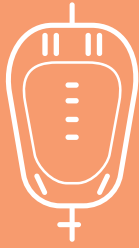

Drainage

Each medical device in our urine drainage range has been designed with a range of features to meet individual patient and clinical requirements in line with NICE² and EPIC¹ best practice guidelines.

IN THIS RANGE

Ugo Foley Catheter Kit	12
Ugo Leg Bag	14
Ugo Pro Leg Bags	16
Ugo 2L Drainage Bags	18
Ugo Sheath	20
Ugo Catheter Valve	22
Ugo Stand with Dignity	24



Our **Ugo Foley Catheter** is a 100% silicone 2-way **Foley catheter**, licensed for transurethral and suprapubic catheterisation for up to 90 days. Each **Ugo Foley Catheter Kit** contains a **Ugo Foley Catheter** and is supplied with a **Ugo 10ml Syringe** for deflation of the Foley catheter balloon and a 10ml **OptiPure** prefilled sterile water syringe for inflation of the Foley catheter balloon.

The **Ugo Foley Catheter** features a universal stepped connector, giving healthcare professionals the freedom to select whichever leg bag or catheter valve they've assessed as being most suited for their patient.

Ugo Foley Catheter Kit

PR CODE	PR NAME	SIZE	OUTER CASE
UFC12	Ugo Foley Catheter Kit	12CH	Case of 12
UFC14	Ugo Foley Catheter Kit	14CH	Case of 12
UFC16	Ugo Foley Catheter Kit	16CH	Case of 12
UFC18	Ugo Foley Catheter Kit	18CH	Case of 12

UFC12 features a non-integrated balloon.

Key features:

- Supplied in a kit with an empty syringe for deflation and sterile water syringe for inflation of the Foley catheter balloon.
- Suitable for transurethral and suprapubic catheterisation.
- Made from 100% silicone.
- Latex-free.
- Integrated 5-10ml balloon on catheters 14CH, 16CH and 18CH.
- Standard length (420mm) suitable for adult male and female catheterisation.
- Suitable for use for up to 90 days.
- X-ray detectable line.
- Includes 3 catheter record stickers to support best practice.





Our **Ugo Leg Bags** are designed for attaching to suprapubic or transurethral catheters or male urinary sheaths, allowing continued use for up to 7 days.

All **Ugo Leg Bags** are supplied sterile, individually pouch packaged when in boxes of 10. Featuring a needle-free sample port for convenient collection of urine, non-return valve to prevent urine going back up the drainage tubing and soft kink-free silicone drainage tubing to create a secure attachment to the 2 litre drainage bag (link drainage) and to reduce the risk of kinking while you sleep.

Every feature of the **Ugo Leg Bag** has been carefully developed with the wearer in mind. Even down to the soft PVC used for the leg bag itself, specifically chosen to reduce noise whilst filling and moving. Our **Ugo Leg Bags** are available in a wide range of formats.

Ugo Leg Bags (Community) - Fabric backing – 1 pair of **Ugo Fix Leg Bag Straps** in every box

PR CODE	PR NAME	BAG VOL	TUBE LENGTH	TAP DESIGN	INNER PACK	OUTER CASE
1C/500/ST/L	Ugo 1C	500ml	Short Tube	Lever Tap	Box of 10	Case of 80
2C/500/LT/L	Ugo 2C	500ml	Long Tube	Lever Tap	Box of 10	Case of 80
3C/500/ST/T	Ugo 3C	500ml	Short Tube	T Tap	Box of 10	Case of 80
4C/500/LT/T	Ugo 4C	500ml	Long Tube	T Tap	Box of 10	Case of 80
5C/350/DI/L	Ugo 5C	350ml	Direct Inlet	Lever Tap	Box of 10	Case of 80
6C/350/ST/L	Ugo 6C	350ml	Short Tube	Lever Tap	Box of 10	Case of 80
7C/750/ST/L	Ugo 7C	750ml	Short Tube	Lever Tap	Box of 10	Case of 80
8C/750/LT/L	Ugo 8C	750ml	Long Tube	Lever Tap	Box of 10	Case of 80

Ugo Leg Bags (Hospital) - Non fabric backing – 1 pair of **Ugo Fix Leg Bag Straps** in every pouch

PR CODE	PR NAME	BAG VOL	TUBE LENGTH	TAP DESIGN	INNER PACK	OUTER CASE
1H/500/ST/L	Ugo 1H	500ml	Short Tube	Lever Tap	Box of 10	Case of 80
2H/500/LT/L	Ugo 2H	500ml	Long Tube	Lever Tap	Box of 10	Case of 80
3H/500/ST/T	Ugo 3H	500ml	Short Tube	T Tap	Box of 10	Case of 80
4H/500/LT/T	Ugo 4H	500ml	Long Tube	T Tap	Box of 10	Case of 80



Key features:

- Suitable for use with a suprapubic or transurethral catheter or male urinary sheath.
- Transparent universal stepped connector for secure fitting into all catheters or sheaths.
- Needle-free sample port, which locks the syringe for safe sample collection.
- Non-return valve to prevent urine going back up the drainage pipe.
- Reinforced eyelets to ensure level positioning when fixated using the **Ugo Fix Leg Bag Straps**.
- Natural leg shape design.
- Durable but soft medical-grade plastic material for strength and reduced noise when moving.



Lever Tap



T Tap

Designed to be used with either an indwelling urethral or suprapubic catheter.



Our **Ugo Pro Leg Bags** are designed for attaching to suprapubic or transurethral catheters or male urinary sheaths, allowing continued use for up to 7 days.

All **Ugo Pro Leg Bags** are supplied sterile, individually pouch packaged with a pair of disposable nitrile gloves and an apron, in boxes of 10.

Featuring a needle-free sample port for convenient collection of urine by healthcare professionals, non-return valve to prevent urine going back up the drainage tubing and soft kink-free silicone drainage tubing to create a secure attachment to the 2 litre drainage bag (link drainage) and to reduce the risk of kinking whilst asleep.

Note: Disposable gloves are not an alternative to hand washing. Ensure you continue to follow hand hygiene guidance before and after changing your patients' **Ugo Pro Leg Bag**.

Ugo Pro Leg Bags .

PR CODE	PR NAME	BAG VOL	TUBE LENGTH	TAP DESIGN	INNER PACK	OUTER CASE
3410	Ugo Pro 1C	500ml	Short Tube	Lever Tap	Box of 10	Case of 80
3411	Ugo Pro 2C	500ml	Long Tube	Lever Tap	Box of 10	Case of 80
3412	Ugo Pro 3C	500ml	Short Tube	T Tap	Box of 10	Case of 80
3413	Ugo Pro 4C	500ml	Long Tube	T Tap	Box of 10	Case of 80



Key features:

- Each **Ugo Pro Leg Bag** is supplied sterile, pouch packaged, with a pair of nitrile gloves and a disposable apron.
- Suitable for use with a suprapubic or transurethral catheter or male urinary sheath.
- Transparent universal stepped connector for secure fitting into all catheters or sheaths.
- Needle-free sample port, which locks the syringe for safe sample collection.
- Non-return valve to prevent urine going back up the drainage pipe.



- Reinforced eyelets to ensure level positioning when fixated using the **Ugo Fix Leg Bag Straps**.
- Natural leg shape design.
- Durable but soft medical-grade plastic material for strength and reduced noise when moving.



Lever Tap



T Tap



Our **Ugo 2L Drainage Bags** are available in sterile and non-sterile variations. Both sterile and non sterile bags are suitable for connecting to a leg bag, as part of a link drainage system. Only our sterile **Ugo 2L Drainage Bags** are suitable for attaching directly to suprapubic or transurethral catheters or male urinary sheaths.

All **Ugo 2L Drainage Bags** are single-use. However, our **Ugo 12** and **Ugo 13** sterile drainage bags can remain in place for up to 7 days.

Ugo 2L Drainage Bags

PR CODE	PR NAME	DESCRIPTION	CLINICAL USE	INNER PACK	OUTER CASE
9/NS/SU/DO	Ugo 9	Non drainable (tear corner to drain) Non sterile 90cm inlet tubing.	Single-use link drainage only	Bag of 10	Case of 80
10/NS/SU/DO	Ugo 10	Single-use T Tap (can only drain once) Non sterile 90cm inlet tubing.	Single-use link drainage only	Bag of 10	Case of 80
11/S/SU/DO	Ugo 11	Single-use T Tap (can only drain once) Sterile pouch packaged. 90cm inlet tubing. 90cm inlet tubing.	Attaches directly to catheter or single-use link drainage	Box of 10	Case of 80
12/S/D/L	Ugo 12	Lever Tap (continuous drainage) Sterile pouch packaged. 90cm inlet tubing.	Attaches directly to catheter for up to 7 days with needle free sample port	Box of 10	Case of 80
13/S/D/T	Ugo 13	T Tap (continuous drainage) Sterile pouch packaged. 90cm inlet tubing.	Attaches directly to catheter for up to 7 days with needle free sample port	Box of 10	Case of 80

Ugo 2L Xtra Drainage Bags - Drainage bags with a longer inlet tube.

PR CODE	PR NAME	DESCRIPTION	CLINICAL USE	INNER PACK	OUTER CASE
3300	Ugo 10xtra	Single Use T Tap (can only drain once) - Non sterile. 120cm inlet tubing.	Single use link drainage only	Bag of 28	Case of 112



Key features:

- Soft, kink-resistant inlet tubing to ensure leak-free continuous drainage.
- Non-return valve to prevent urine going back up the drainage pipe.
- Reinforced eyelets to ensure level hanging on the bag stand.
- Needle-free sprung sample port, which locks the syringe for safe sample collection (Ugo 12 and Ugo 13).
- Transparent universal stepped connector for secure fitting into all indwelling catheters or leg bag silicone outlet tubes for overnight link drainage.
- Provided in a tamper-resistant, easy to store, resealable bag (Ugo 9, Ugo 10 and Ugo 10 Xtra).

Designed to be used for link drainage and/or continuous free drainage.



The **Ugo Sheath** is available in a standard format and a pop-on format, each available in five size options.

Compatible with all catheter leg bags, to provide a leak-free connection, so you can go about your day with confidence.

The **Ugo Sheath** has been created using minimal plastic and with cardboard packaging that can be fully recycled.

Ugo Sheath

PR CODE	DESCRIPTION	INNER PACK	OUTER CASE
3100	Ugo Sheath Standard - 25mm	Box of 28	Case of 112
3101	Ugo Sheath Standard - 28mm	Box of 28	Case of 112
3102	Ugo Sheath Standard - 30mm	Box of 28	Case of 112
3103	Ugo Sheath Standard - 32mm	Box of 28	Case of 112
3104	Ugo Sheath Standard - 36mm	Box of 28	Case of 112
3105	Ugo Sheath Standard - 40mm	Box of 28	Case of 112
3106	Ugo Sheath Pop-on - 25mm	Box of 28	Case of 112
3107	Ugo Sheath Pop-on - 28mm	Box of 28	Case of 112
3108	Ugo Sheath Pop-on - 30mm	Box of 28	Case of 112
3109	Ugo Sheath Pop-on - 32mm	Box of 28	Case of 112
3110	Ugo Sheath Pop-on - 36mm	Box of 28	Case of 112
3111	Ugo Sheath Pop-on - 40mm	Box of 28	Case of 112



Key features:

- Made of silicone.
- Latex-free.
- Easy to apply.
- Discreetly packaged.
- Each box is fully recyclable.
- Kinder to the environment. Made with less plastic.
- Pubic hair protector included for more comfortable application.
- Resealable travel pouch provided with every box.





Easy to use, discreet and comfortable to wear, and designed in detail for safety and security. The **Ugo Catheter Valve** gives patients the confidence they need to carry on living independent lives. Using a catheter valve as an alternative to urine drainage bags may improve experience and outcomes for suitable patients.

The **Ugo Catheter Valve** is a small, sterile, single-use drainage valve which is used to allow intermittent urinary drainage from an indwelling urethral or suprapubic catheter.

The **Ugo Catheter Valve** features a universal stepped connector which is connected directly to the catheter, and a lever for controlled manual urine drainage. A **Ugo Catheter Valve** can remain in place for up to 7 days.

*"I use the **Ugo Catheter Valve** with my suprapubic catheter. It grips extremely well, and doesn't slip like some other valves which I have tried."*


Christine P, London, UK

Ugo Catheter Valve

PR CODE	DESCRIPTION	INNER PACK	OUTER CASE
3000	Ugo Catheter Valve	Box of 5	Case of 200

Key features:

- Offers a discreet alternative to urinary drainage bags for people assessed and deemed suitable.
- Allows the bladder to fill and hold urine.
- Can remain in place for up to 7 days.
- Suitable for connecting to a drainage bag for link drainage if required.
- Safety lever feature to prevent accidental opening.
- Universal stepped connector for secure fit
- Lever handle designed to allow easy single-handed operation.
- Contrasting colour safety lever for easy visual identification.

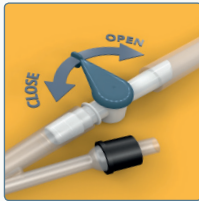



Thoughtful and creative, the **Ugo Stand with Dignity** is provided with a dignity cover to discreetly cover the drainage bag if required. The **Ugo Stand with Dignity** is designed for convenience and discretion, with multiple positioning options. It can be positioned free standing on the floor, hooked onto a bed frame or wheelchair, or slotted between the mattress and divan.

The **Ugo Stand with Dignity** can also be easily moved if required using the handle. The **Ugo Stand with Dignity** is reusable and can be easily wiped down or disinfected in a sleuth.

*"I'm so grateful that **Optimum Medical** have created the **Ugo Stand with Dignity**. The removable dignity cover features a 'viewing strip' so you can see how full your bag is getting without having to have anything personal on show. Thank you for giving me back my privacy and dignity – this is so important to me."*

Margaret Balch, Birmingham, UK

Ugo Stand with Dignity

PR CODE	DESCRIPTION	INNER PACK	OUTER PACK
3024	Ugo Stand with Dignity	Box of 1	Case of 4

Key features:

- Dignity cover provides discreet urine drainage.
- Purpose designed to reduce the risk of infection by ensuring the tap doesn't touch the ground.
- Variety of usage options including free standing or hooked onto a bed frame or wheelchair.





Fixation Devices

Purpose designed fixation devices for safety and effectiveness to reduce the risk of trauma, pain and infection.

NICE² guidelines state that 'best practice suggests catheters are secured to avoid trauma'.

IN THIS RANGE

Ugo Fix Gentle (catheter clip)	28
Ugo Fix Catheter Strap	30
Ugo Fix Sleeve (leg bag holder)	32
Ugo Fix Bag Belt	34
Multi Fix Gentle	35



Our **Ugo Fix Gentle** is a non-sterile, medical device, designed to comfortably adhere to the skin, providing fixation for a two-way indwelling urinary catheter.

The **Ugo Fix Gentle** features a rotating clip which is designed to hold the catheter away from the skin as it rotates. The catheter can move naturally with the motion of the body, reducing the risk of pain or discomfort within the urethra or at the bladder neck from the catheter moving or pulling inadvertently.

There is no skin preparation required prior to applying the **Ugo Fix Gentle**, and it can be gently removed without the use of solutions or chemicals. The **Ugo Fix Gentle** catheter clip features silicone binding, which can hold in place suprapubic or transurethral catheters, ranging in size from 8CH to 20CH.

The benefits of our silicone re-adhesive gel pad:


- Provides immediate, gentle adhesion to the skin.
- Easy to apply and easy to remove.
- No skin preparation required before application.
- No need for solutions to dissolve the adhesive when removed.
- Leaves no sticky residue.
- Showerproof, breathable and translucent.
- Suitable for up to 7 days continuous use depending on users lifestyle, activity levels and ability to self-care.

Ugo Fix Gentle (catheter clip)

PR CODE	DESCRIPTION	INNER PACK	OUTER CASE
3004	Ugo Fix Gentle (catheter clip)	Box of 5	Case of 100

The benefits of our rotating catheter clip:

- Internal silicone bonding to effectively grip catheters ranging in size from 8CH to 20CH.
- Revolving clip allows natural movement of the catheter with the body.
- Raised arrow to indicate direction of catheter insertion into the clip.
- Easy-open catch.
- Positive sensory click to indicate closure.
- Holds tubing away from skin to reduce discomfort.
- Contrasting colour clip for easy visual identification.

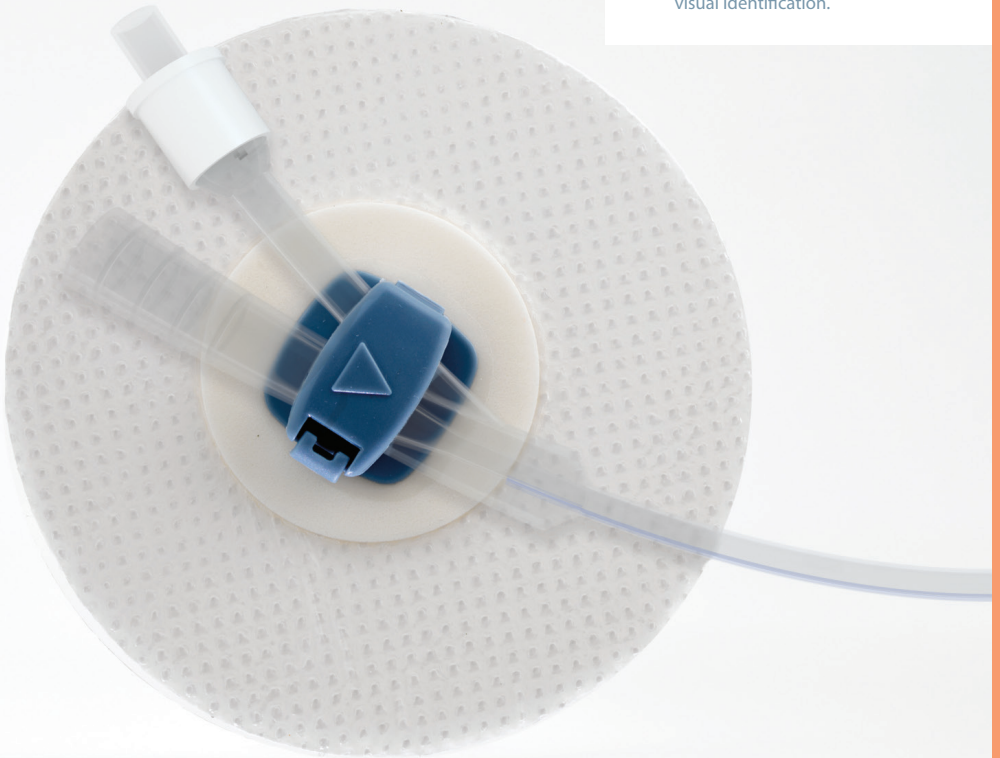

With silicone technology for atraumatic patient experience

ugofix™

catheter strap



Secure, reliable and simple to use. The **Ugo Fix Catheter Strap** provides comfortable support of the indwelling urinary catheter, giving your patients the freedom to move around freely and safely.

Available in five sizes, all of which can be cut to size.

Suitable for use on the abdomen or thigh.



To assist with patient's daily wash and wear cycle, we provide a laundry bag in each box of **Ugo Fix Catheter Straps**.

*"The **Ugo Fix Catheter Strap** made life more comfortable for my husband while he had a catheter. The previous straps he used cut into his leg whereas the **Ugo Fix Catheter Strap** was wider and longer, so we could easily cut to size, and we found the fastening efficient and secure. They washed well and dried quickly".*

Nicola Bradbury, Reading, UK

Ugo Fix Catheter Strap

PR CODE	DESCRIPTION	INNER PACK	OUTER CASE
3001	Ugo Fix Catheter Strap Extra short – up to 35cm	Box of 5	Case of 50
3002	Ugo Fix Catheter Strap Short – up to 45cm	Box of 5	Case of 50
3003	Ugo Fix Catheter Strap Medium – up to 80cm	Box of 5	Case of 50
3009	Ugo Fix Catheter Strap Long – up to 130cm	Box of 5	Case of 50
3010	Ugo Fix Catheter Strap Extra Long – up to 180cm	Box of 5	Case of 50



Key features:

- Designed to secure the catheter to prevent trauma within the urethra or bladder neck.
- All five lengths are fully adjustable and can be cut to fit.
- Silicone technology for comfortable, yet secure grip to skin.
- Reusable and washable (provided with laundry bag).
- Ideally complemented by the **Ugo Fix Sleeve** and **Ugo Fix Gentle**.
- Contrasting tab for easy visual identification.

Designed to support urethral catheters, suprapubic catheters and sheath systems. Fixes on either the catheter or leg bag tubing.



Durable, comfortable and extra supportive. The **Ugo Fix Sleeve** is designed to support leg bags attached to indwelling catheters and sheaths.

The **Ugo Fix Sleeve** is available in five sizes, each of which can support a leg bag from 350ml to 750ml in volume.

The **Ugo Fix Sleeve** is made from elasticated fabric for added user comfort, and features external seams and a reinforced top and bottom band.



To assist with patient's daily wash and wear cycle, we provide a laundry bag in each box of **Ugo Fix Sleeves**.

*"The **Ugo Fix Sleeve** leg bag holder has given my partner enough confidence to wear shorts this summer - something he usually does but was feeling he might not be able to do with a leg bag. The **Ugo Fix Sleeve** has made this option available to him and given him the greatest confidence boost!"*

Martin Otter, Devon.

Ugo Fix Sleeve (leg bag holder)

PR CODE	DESCRIPTION	INNER PACK	OUTER CASE
3005	Ugo Fix Sleeve Small – 24cm-39cm	Box of 4	Case of 80
3006	Ugo Fix Sleeve Medium – 36cm-55cm	Box of 4	Case of 80
3007	Ugo Fix Sleeve Large – 40cm-70cm	Box of 4	Case of 80
3008	Ugo Fix Sleeve Extra Large – 65cm-90cm	Box of 4	Case of 80
3027	Ugo Fix Sleeve XXL - 75cm-105cm	Box of 4	Case of 80



Designed for supporting leg bags attached to indwelling catheters and sheaths

Key features:

- Suitable for supporting leg bags attached to indwelling catheters and sheaths.
- Available in five size options, each of which can support a leg bag volume of 350ml to 750ml.
- No need to disconnect leg bag from catheter when changing sleeve.
- External seam to protect skin and increase comfort.
- Colour coded for easy reordering.
- Reusable and washable (provided with a laundry bag).
- Ideally complemented by the **Ugo Fix Gentle** and **Ugo Fix Catheter Strap**.

ugofix™

bag belt



The **Ugo Fix Bag Belt** is a non-sterile medical device designed for supporting a urinary leg bag that is attached to an indwelling urinary catheter or male urinary sheath.

The **Ugo Fix Bag Belt** is 70cm long and is suitable for users with a leg circumference up to 65cm. It is adjustable and can be cut down to fit if needed.

The **Ugo Fix Bag Belt** is washable and reusable and can support urinary leg bags from 350ml to 750ml volume

Key features:

- Designed for supporting leg bags ranging from 350ml to 750ml.
- Wider and padded strap for additional comfort when in use.
- Holds leg bag away from skin for additional comfort.
- One size (70cm in length), can be cut to size.
- Washable and reusable.
- Contrasting coloured tabs for easy visual identification.



Ugo Fix Bag Belt

PR CODE	DESCRIPTION	INNER PACK	OUTER CASE
3012	Ugo Fix Bag Belt	Box of 3	Case of 60

multifix™

gentle



The **Multi Fix Gentle** is a non-sterile medical device designed to secure indwelling urinary catheters and medical tubes up to 10mm in diameter. The **Multi Fix Gentle** is intended to secure the catheter or medical tubing, providing comfort for the user and reducing the risk of trauma as a result of inadvertent movement of the catheter or medical tubing.

The **Multi Fix Gentle** features a tubing securement tab which is designed to hold the catheter or tubing in place without causing damage to the tubing. The tubing securement tab is mounted on a re-adhesive silicone gel strip, which is applied directly to the skin without the need for any skin preparation.

The **Multi Fix Gentle** is available in two sizes, small (110mm x 25mm) and large (140mm x 40mm).





Key features:

- Secures urinary catheters and medical tubing up to 10mm in diameter.
- Suitable for securing medical tubing and urinary catheters.
- Easy to apply and easy to remove.
- No skin preparation required.
- Leaves behind no sticky residue.
- For up to 7 days continuous use.
- Available in two sizes.

Multi Fix Gentle - universal tubing securement device

PR CODE	DESCRIPTION	INNER PACK	OUTER CASE
3030	Multi Fix Gentle (large)	Box of 10	Case of 200
3031	Multi Fix Gentle (small)	Box of 10	Case of 200

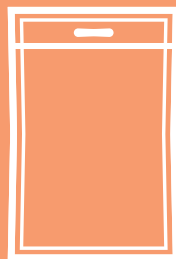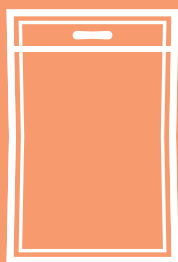

Catheter drainage and management packs

Smart management systems created to reduce product shortages and improve efficiency.

IN THIS RANGE

Ugo 4 Weeks	38
Ugo Home	40



The **Ugo 4 Weeks** is a 4 week supply of essential catheter drainage and fixation medical devices, all contained in a compact, discreet box.

One Box

With protective and lightweight packaging, the **Ugo 4 Weeks** is easy to store in a person's home or care home setting, whilst also reducing stockpiling and wastage.

One Code

The **Ugo 4 Weeks** makes re-ordering simple, with just one code required.

One Prescription

It's easy and efficient to prescribe, saving time for everyone involved.

One Month

The **Ugo 4 Weeks** effectively prevents product shortages and improves efficiency, improving continuity of care.

Ugo 4 Weeks

PR CODE	PRODUCT NAME	LEG BAG DESCRIPTION
3014	Ugo 4 Weeks 1C	(500ml, Short Tube, Lever Tap)
3015	Ugo 4 Weeks 2C	(500ml, Long Tube, Lever Tap)
3016	Ugo 4 Weeks 3C	(500ml, Short Tube, T Tap)
3017	Ugo 4 Weeks 4C	(500ml, Long Tube, T Tap)



Each box contains:

- 4 x fabric backed **Ugo Leg Bags** (1C/2C/3C/4C).
- 28 x **Ugo 10 2L Bags** with single use T Tap.
- 1 x **Ugo Hanger**.
- 4 x pairs of **Ugo Fix Leg Bag Straps**.
- 4 x **Ugo Fix Gentle** catheter clips.





Convenient and compact catheter drainage and management system, containing one week's supply of products. Often referred to as our "Discharge Pack", the **Ugo Home** is ideal for use when patients are discharged from hospital with an indwelling catheter. Each product inside the **Ugo Home** is available separately.

Each pack contains:

- 1 x fabric backed **Ugo Leg Bag** (**Ugo 1C / Ugo 2C / Ugo 3C / Ugo 4C**).
- 7 x **Ugo 10 2L Bags** with single use T Tap.
- 1 x **Ugo Hanger**.
- 1 x pair of **Ugo Fix Leg Bag Straps**.
- 1 x **Ugo Fix Gentle** (catheter clip).

Ugo Home - 7-day discharge pack

PR CODE	DESCRIPTION	INNER PACK	OUTER CASE
3018	Ugo Home 1C	1 Bag	8 Bags
3019	Ugo Home 2C	1 Bag	8 Bags
3020	Ugo Home 3C	1 Bag	8 Bags
3021	Ugo Home 4C	1 Bag	8 Bags



x1



x7



x1



x1 pair



x1



x1





Catheterisation Jelly

A range of catheterisation jelly, supporting patient-centred care by offering a choice.

NICE² guidelines state that 'urethral trauma and discomfort will be minimised by using an appropriate sterile, single-use lubricant or anaesthetic gel'.

IN THIS RANGE

OptiLube Active	45
OptiLube CHG Free	46
OptiLube syringes	47
Choosing the right catheterisation jelly for your patient	48

Our catheterisation jelly range

We believe in choice. Which is why we offer a range of catheterisation jelly pre-filled syringes, designed specifically for use during catheterisation procedures.

All our catheterisation jelly syringes are supplied sterile inside and out and designed to fall from their packaging onto a sterile field, supporting an aseptic technique.

- Each available in 6ml (female) and 11ml (male).
- Manufacturing date, expiry date and LOT number are printed on the tear off strips.
- Provided in a familiar and recommended syringe format.
- Easy-to-prime syringe for controlled application.
- Rounded tip for comfort when inserted into the urethra.

For more information on any product in our catheterisation jelly range, email enquiries@optimummedical.co.uk or call +44 (0) 845 643 5479

optilube.TM active



6ml

11ml

Containing local anaesthetic (Lidocaine 2%) and antiseptic (Chlorhexidine Gluconate 0.05%).



OptiLube Active

PR CODE	PR NAME	SIZE	INNER PACK	OUTER PACK
1160	OptiLube Active	6ml	10	200
1161	OptiLube Active	11ml	10	200

optilube.™

active CHG free



6ml

11ml

Containing local anaesthetic (Lidocaine 2%). Ideally suited for patients with a contraindication or sensitivity to Chlorhexidine Gluconate (CHG).



Optilube Active CHG Free

PR CODE	PR NAME	SIZE	INNER PACK	OUTER PACK
1167	Optilube Active CHG Free	6ml	10	200
1168	Optilube Active CHG Free	11ml	10	200

optilube.™



6ml

11ml

Containing no active ingredients. Ideally suited for patients with a contraindication or sensitivity to Lidocaine or Chlorhexidine Gluconate (CHG).



OptiLube

PR CODE	PR NAME	SIZE	INNER PACK	OUTER PACK
1125	OptiLube Syringe	6ml	25	250
1126	OptiLube Syringe	11ml	25	250

Choosing the right catheterisation jelly for your patient

To making an informed choice for your patient, you need to know what options are available to you. We offer three choices of catheterisation jelly, each available in 6ml (female) or 11ml (male) pre-filled syringes.

This is a visual guide only. The choice of which catheterisation jelly to use should be made following an individual patient assessment and in line with local healthcare policies and procedures.

For more guidance on selecting the right catheterisation jelly from our range, email enquiries@optimummedical.co.uk or call +44 (0) 845 643 5479

Suprapubic

We recommend 6ml **OptiLube** for all routine suprapubic catheter changes because of the swift change of catheter required during this procedure.*

Urethral

Does the patient have any sensitivities or allergies to Chlorhexidine or are there any contraindications to its use?

YES

NO

Does the patient have any sensitivities or allergies to Lidocaine or are there any contraindications to its use?

Does the patient have any sensitivities or allergies to Lidocaine or are there any contraindications to its use?

YES

NO

OptiLube is the ideal choice for the patient. Available in 11ml (male) or 6ml (female) formats.

OptiLube is the ideal choice for the patient. Available in 11ml (male) or 6ml (female) formats.

YES

NO

You can use **OptiLube** or **OptiLube Active CHG Free**. Both are available in 11ml (male) or 6ml (female) formats.

You can use **OptiLube**, **OptiLube Active** or **OptiLube Active CHG Free**. All are available in 11ml (male) or 6ml (female) formats.

Sustainability

Over the years we've noticed a lack of attention in the healthcare industry towards sustainability. We felt we could do something about that – with even the smallest changes growing into great things. Change starts with our people, so we're making sure that everyone at **Optimum Medical** understands the role they play in making things better for our environment and our communities.

We've started to lay the groundworks for a greener future at our offices and warehouse, with electric vehicle charging points, switching to LED lighting and purchasing our energy through renewable suppliers. Our packaging is procured responsibly, using FSC approved sources and purchasing recycled materials wherever we are able.

And that's just the beginning.

We will continue to strive to make adjustments and improvements in every area of our business. For more information on the steps we are taking, please get in touch with us on +44 (0) 845 643 5479 or email sustainability@optimummedical.co.uk



NOTES
